
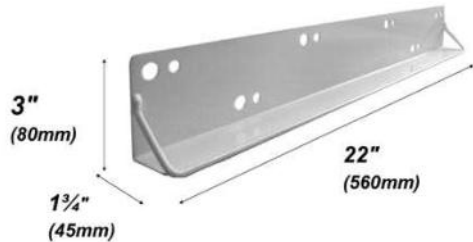


Myco Wall Bracket (MWB)

Allows a Swivel Bather to be used on baths with a narrow ledge on the wall side.

The Myco Wall Bracket is designed to support a Swivel Bather when the lip on the wall side of a bath is less than the recommended 1 1/4" (30mm).

The Stainless Steel construction allows it to be used for all Myco Swivel Bather models. It can also be used for other products requiring a ledge on the wall side of the bath (e.g. Bath Boards)

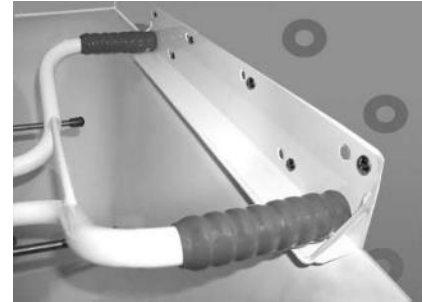

The Wall Bracket is located so that it rests on the edge of the bath as shown and is secured to the wall using suitable fixings (not supplied).